


English

OPERATING MANUAL

Water Baths

TW2
TW8
TW12
TW20

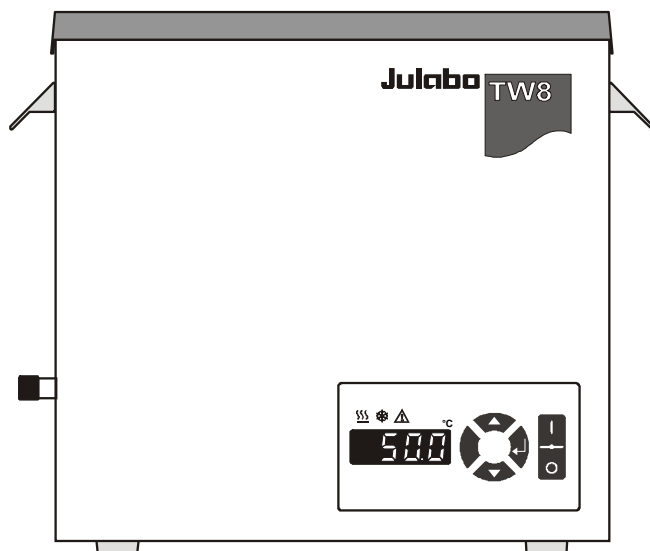

Julabo
THE TEMPERATURE CONTROL COMPANY

JULABO Labortechnik GmbH
77960 Seelbach / Germany

☎ +49 (0) 7823 / 51-0

📄 +49 (0) 7823 / 24 91

✉ info@julabo.de

🌐 www.julabo.de

1.951.6062BE4 03/10

19516062.doc

Version: 12.04.10

Congratulations!

You have made an excellent choice.

JULABO thanks you for the trust you have placed in us.

This operating manual has been designed to help you gain an understanding of the principles of operating and possibilities of our water baths. For optimum utilization of all functions, we recommend that you thoroughly study this manual prior to beginning operation.

Quality Management System



The JULABO **Quality Management System**:

Development, production and distribution of temperature application instruments for research and industries conform to the requirements according to DIN EN ISO 9001:2000.

Certificate Registration No. 01 100044846

Unpacking and checking

Unpack the water baths and accessories and check for damages incurred during transit.

These should be reported to the responsible carrier, railway, or postal authority, and a request for a damage report should be made. These instructions must be followed fully for us to guarantee our full support of your claim for protecting against loss from concealed damage.

The form required for filing such a claim will be provided by the carrier.

Printed in Germany

Changes without prior notification reserved

TABLE OF CONTENTS

Operating manual.....	4
1. Intended use	4
1.1. Description.....	4
2. Operator responsibility – Safety recommendations.....	4
2.1. Disposal.....	6
2.2. Use in Sous-Vide cooking.....	7
2.3. EC Conformity	8
2.4. Warranty conditions	8
2.5. Technical specifications	9
Operating instructions	12
3. Operating controls and functional elements.....	12
4. Safety notes for the user.....	13
4.1. Explanation of safety notes.....	13
4.2. Explanation of other notes	13
5. Safety recommendations	13
6. Preparations.....	15
6.1. Installation.....	15
6.2. Bath fluids	16
6.3. Filling / Draining	17
6.4. Maintaining a constant water level / Counter cooling.....	18
6.5. Accessoires	19
7. Operating procedures	20
7.1. Power connection	20
7.2. Switching on / Start - Stop	21
7.3. Setting the temperatures.....	21
8. Troubleshooting guide / Error messages	22
9. Cleaning / repairing the unit	22

Operating manual

1. Intended use

JULABO water baths have been designed for temperature application to specific fluids in a bath tank.




JULABO water baths are not conceived for direct temperature application to food and luxury articles or pharmaceutical and medico-technical products. Direct temperature application means: Unprotected contact of the object with the bath medium (bath fluid).

1.1. Description



☑ The water baths are operated via the splash-proof keypad. The implemented microprocessor technology allows to set and to store the setpoint that can be indicated on the LED temperature display.

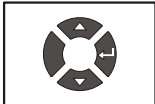

☑ The PID temperature control adapts the heat supplied to the thermal requirements of the bath.



☑ The water baths conforms to the relevant requirements specified by European guidelines.

☑ The models TW8, TW12, TW20 provide a drain screw and handles for easy relocation.

2. Operator responsibility – Safety recommendations

The products of JULABO Labortechnik GmbH warrant a safe operation if installation, operation and maintenance is carried out according to common safety regulations. This section informs you about potential dangers that may arise from operating the circulator and also mentions the most important safety precautions





The operator is responsible for the qualification of the personnel operating the units.

- The operator should be constantly informed about the dangers involved with their job activities as well as preventive actions.
- Make sure all persons expected to carry out operation, installation and maintenance of the unit read and understand the safety information and operating instructions.

- When using hazardous materials, the circulator may only be operated by persons that are absolutely familiar with these materials and the unit. These persons must be fully aware of possible risks.

If you have any questions concerning the operation of your unit or the information in this manual, please contact us!

Contact JULABO Labortechnik GmbH
Eisenbahnstraße 45
77960 Seelbach / Germany

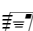 +49 (0) 7823 / 51-0
 +49 (0) 7823 / 24 91
 info@julabo.de
 www.julabo.de

Safety recommendations for the operator

- You received a product conceived for industrial use. Nevertheless, avoid strikes to the housing, vibrations, damages to the keypad foil (keys, display) or contamination.
- Make sure the product is regularly checked for proper condition. Regularly check (at least every 2 years) the proper condition of the mandatory, warning, prohibition and safety labels.
- Take care that the mains supply features a low impedance to avoid any negative affects on the instrument being operated in the same mains.
- This unit is designed for operation in a controlled electromagnetic environment. This means that transmitting devices (e.g. cellular phones) should not be used in the immediate vicinity. Magnetic radiation may influence other units with components susceptible to magnetic fields (e.g. a monitor). We recommend to keep a minimum distance of 1 m.
- Permissible ambient temperature: max. 40 °C, min. 5 °C.
- Permissible relative air humidity: 50 % (40 °C).
- Do not store in an aggressive atmosphere. Protect from contaminations.
- Do not expose to sunlight.

Appropriate Operation

Only qualified personnel is authorized to perform configuration, installation, maintenance and repairs of the water bath.

Routine operation can also be carried out by untrained personnel who should however be instructed by trained personnel.

Use:

Insufficient ventilation may result in the formation of explosive mixtures. Only use the unit in well ventilated areas. The unit is not for use in explosive atmosphere.




JULABO shaking water baths have been designed for temperature application to water in a bath tank.

The bath may **not** be filled with flammable materials. Fire hazard!


Only use non-acid and non corroding bath fluids.

Operator responsibility – Safety recommendations


When using hazardous materials or materials that could become hazardous, **the operator must** affix the enclosed safety labels **(1 + 2)** to the front of the unit so they are highly visible:

1		Warning label W00: Colors: yellow, black Danger area. Attention! Observe instructions. (operating manual, safety data sheet)
2		Mandatory label M018: Colors: blue, white Carefully read the user information prior to beginning operation. Scope: EU
or		
2		Semi S1-0701 Table A1-2 #9 Carefully read the user information prior to beginning operation. Scope: USA, NAFTA

Particular care and attention is necessary because of the wide operating range. There are thermal dangers: Burn, scald, hot steam, hot parts and surfaces that can be touched.

	Warning label W26: Colours: yellow, black Hot surface warning. (The label is put on by JULABO)
--	--

2.1. Disposal

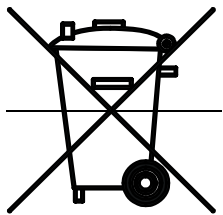

Valid in EU countries

Directive 2002/96/EC of the European Parliament and of the Council of 27 January 2003 on waste electrical and electronic equipment (WEEE).

This directive requires electrical and electronic equipment marked with a crossed-out trash can to be disposed of separately in an environmentally friendly manner.

Contact an authorized waste management company in your country.

Disposal with household waste (unsorted waste) or similar collections of municipal waste is not permitted!

2.2. Use in Sous-Vide cooking

Some JULABO product may be used for low temperature cooking.

Sous-Vide cooking should only be used by trained professionals.

For proper food safety it is important to maintain precise and consistent conditions and a hygienic environment.

It is vitally important to maintain procedures that minimize the risk from anaerobic bacteria, e.g. *costridium botulinum*.

Please make sure to follow the Sous-Vide techniques of your trained professional or cooking publication and we remind you to consider the following:

- Use high quality ingredients from known suppliers;
- Separate raw and cooked food areas from each other;
- Establish minimum cooking temperature for meat and fish;
- Adhere to the maximum storage temps and times suggested by your teacher or literature;
- Label contents with date of preparation as well as expiration date;
- Operate in a sterile environment including, utilization of sterile gloves when filling the Sous-Vide bags;
- Use separate vacuums for raw and cooked materials.

Additionally, it is important that your kitchen staff employ an appropriate HACCP (Hazard s Analysis and Critical Control Point Program) to address the movement of food from delivery to service.

Monitor the quality of water in the Sous-Vide bath and change regularly.

Be observant of potential breaks in the cold chain or unnecessary holding time.

Additionally, please make sure to employ the appropriate time to cool cooked food before storage.

To achieve pasteurization, food cooked in the Sous-Vide method should achieve at least 65 °C (149 °F) at its core. Utilization of a temperature probe or other thermometer may be required.

Last, please follow the Sous-Vide cooking regulations mandated by your local Health Bureau or any applicable food regulations for your jurisdiction.

2.3. EC Conformity



The products described in the operating instructions conform to the requirements of the following European guidelines:

Low voltage regulations with respect to legal harmonization of the member countries concerning electric devices for use within certain voltage limits.

EMC guideline with respect to legal harmonization of the member countries concerning electromagnetic compatibility.

Julabo

JULABO Labortechnik GmbH
Eisenbahnstr. 45
77960 Seelbach / Germany

2.4. Warranty conditions

JULABO Labortechnik GmbH warrants its products against defects in material or in workmanship, when used under appropriate conditions and in accordance with appropriate operating instructions

for a period of ONE YEAR.

Extension of the warranty period – free of charge

2 Years Warranty

1Plus Warranty

Registration free of charge on www.julabo.de

With the '1PLUS warranty' the user receives a free of charge extension to the warranty of up to 24 months, limited to a maximum of 10 000 working hours.

To apply for this extended warranty the user must register the unit on the JULABO web site www.julabo.de, indicating the serial no. The extended warranty will apply from the date of JULABO Labortechnik GmbH's original invoice.

JULABO Labortechnik GmbH reserves the right to decide the validity of any warranty claim. In case of faults arising either due to faulty materials or workmanship, parts will be repaired or replaced free of charge, or a new replacement unit will be supplied.

Any other compensation claims are excluded from this guarantee.

2.5. Technical specifications

Water Bath			TW2	TW8
Working temperature range	°C		20 ... 99,9*	20 ... 99,9*
Temperature stability			±0,2	±0,2
Temperature selection			digital	digital
Temperature indication			LED	LED
Resolution	°C		0.1	0.1
Temperature control			PID1	PID1
Heater wattage (at 230 V)	kW		1,0	2,0
Heater wattage (at 115 V)	kW		1,0	1,0
Bath opening (WxL)	cm		15x13	23x27
Bath depth	cm		11	14
Filling volume	liters		1 ... 2	3 ... 8
Overall dimensions (WxDxH) with Makrolon [®] -cover	cm		17x16x26 17x16x37	29x32x28 29x32x44
Weight	kg		3,5	8,5
Ambient temperature	°C		5 ... 40	5 ... 40
Mains power connection 230 V/50-60 Hz	V/ Hz		190-253 / 50-60	190-253 / 50-60
Current draw (at 230 V)	A		4	8
Mains power connection 115 V/60 Hz	V/ Hz		103-127 / 50-60	103-127 / 50-60
Current draw (at 115 V)	A		9	9

All measurements have been carried out at: rated voltage and frequency
 operating temperature: 70 °C ambient temperature: 20 °C bath fluid: water

Technical changes without prior notification reserved.

Water Bath			TW12	TW20
Working temperature range	°C		20 ... 99,9*	20 ... 99,9*
Temperature stability			±0,2	±0,2
Temperature selection			digital	digital
Temperature indication			LED	LED
Resolution	°C		0.1	0.1
Temperature control			PID1	PID1
Heater wattage (at 230 V)	kW		2,0	2,0
Heater wattage (at 115 V)	kW		1,0	1,0
Bath opening (WxL)	cm		35x27	50x30
Bath depth	cm		14	18
Filling volume	liters		5 ... 14	8 ... 26
Overall dimensions (WxDxH) with Makrolon [®] -cover	cm		40x32x28 40x32x44	50x35x32 50x35x49
Weight	kg		8,9	14,2
Ambient temperature	°C		5 ... 40	5 ... 40
Mains power connection 230 V/50-60 Hz	V/ Hz		190-253 / 50-60	190-253 / 50-60
Current draw (at 230 V)	A		8	8
Mains power connection 115 V/60 Hz	V/ Hz		103-127 / 50-60	103-127 / 50-60
Current draw (at 115 V)	A		9	9

All measurements have been carried out at: rated voltage and frequency
operating temperature: 70 °C ambient temperature: 20 °C bath fluid: water

Technical changes without prior notification reserved.

Safety installations according to IEC 61010-2-010:

Excess temperature protection	105 °C - fixed value
Classification according to DIN 12876-1	class I
Alarm message	optical + audible (permanent)

Environmental conditions according to EN 61 010, part 1:

Use only indoor.

Altitude up to 2000 m - normal zero.

Ambient temperature: +5 ... +40 °C (for storage and transportation)

Air humidity:

Max. rel. humidity 80 % for temperatures up to +31 °C,

linear decrease down to 50 % relative humidity at a temperature of +40 °C

Max. mains fluctuation of ± 10 % are permissible.

Protection class according to EN 60 529 IP43

The unit corresponds to Class I

Overvoltage category II

Pollution degree 2



Caution:

The unit is not for use in explosive environment

Standards for interference resistance according to EN 61326-1

This unit is an ISM device classified in Group 1 (using high frequency for internal purposes)


Class A (industrial and commercial range).

Operating instructions





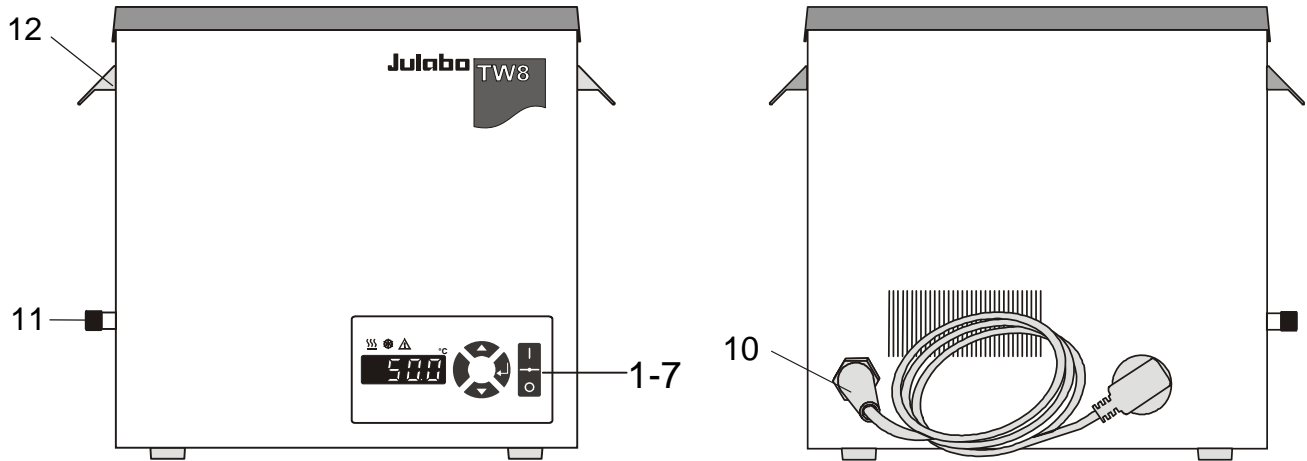


3. Operating controls and functional elements

Front view

Rear view



Example: TW8

- | | | |
|----|---|---|
| 1 |  | <p>Mains power switch, illuminated</p> <p style="text-align: center;"> on
 off </p> |
| 2 |  | <p>Edit keys (increase/decrease setting)</p> |
| 3 |  | <p>Enter key (store)</p> |
| 4 |  | <p>LED temperature display, menu indication</p> |
| 5 |  | <p>Control indicator – Heating</p> |
| 6 |  | <p>Control indicator – Cooling (without function)</p> |
| 7 |  | <p>Control indicator – Alarm</p> |
| 10 | | <p>Mains power cable with plug</p> |
| 11 | | <p>Drainage screw
Connector for liquid level/cooling set (accessory)</p> |
| 12 | | <p>Handle</p> |

4. Safety notes for the user

4.1. Explanation of safety notes



In addition to the safety warnings listed above, warnings are posted throughout the manual. These warnings are designated by an exclamation mark inside an equilateral triangle. „Warning of a dangerous situation (Attention! Please follow the documentation).“
The danger is described according to an alarm keyword.
Read and follow these important instructions.



Warning:

Describes a possibly highly dangerous situation. If this is not avoided, serious injury and danger to life could result.



Caution:

Describes a possibly dangerous situation. If this is not avoided, slight or minor injuries could result.
A warning of possible damage can also be contained in the text.



Notice:

Describes a possibly harmful situation. If this is not avoided, the product or anything in its surroundings can be damaged.

4.2. Explanation of other notes



Note!

Wants to draw your attention to something particular.



Important!

Describes useful information for the operation and the user.

5. Safety recommendations

Follow the safety recommendations to prevent damage to persons or property. Further, the valid safety instructions for working places must be followed.



- ConnOnly connect the unit to a power socket with earthing contact (PE – protective earth)!
- Operation is permitted with **non-flammable** liquids only.
- Place the instrument on an even surface on a pad made of **non-inflammable** material.
- Do not stay in the area below the unit.

Safety recommendations

- Make sure you read and understand all instructions and safety precautions listed in this manual before installing or operating your unit.
- Never operate the unit without bath fluid in the bath.
- Prevent water from penetrating into the hot bath oil.
- Do not drain the bath fluid while it is hot!
Check the temperature of the bath fluid prior to draining (by switching the unit on for a short moment for example).
- Never operate damaged or leaking equipment.
- Always turn off the unit and disconnect the mains cable from the power source before performing any service or maintenance procedures, or before moving the unit.
- Always empty the bath before moving the unit.
- Never operate equipment with damaged mains power cables.
- Condensation that could appear in and on other units near the water bath may result in reduced operating safety.
Be careful when setting up and operating the water bath!



- Some parts of the bath cover may become extremely warm during continuous operation.
When lifting the bath cover, pay attention to hot steam!



- Be careful when touching these parts!
- Use safety glasses!



Use of your JULABO product in Sous-Vide cooking should only be done by someone trained in proper Sous-Vide cooking techniques and safety. For proper food safety it is important to maintain precise consistent conditions and a hygienic environment. Proper procedures minimize the risk from anaerobic bacteria, Julabo makes no warranty and excludes all claims concerning food products prepared under Sous-Vide techniques.


6. Preparations

6.1. Installation



Caution:

The unit is not for use in explosive environment.

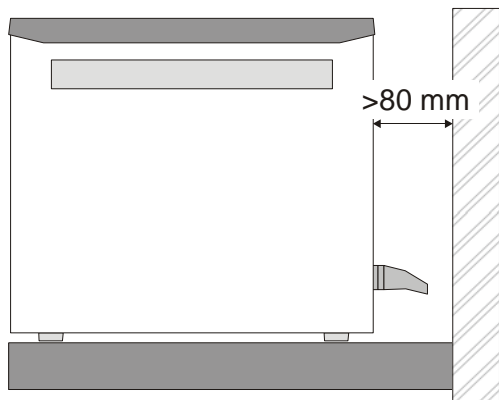

Place the waterbath in an upright position.

Keep a wall distance of minimum 80 mm.



Caution:

Potential hazards from the samples

Proper use of shaking water baths includes immersion of samples contained in test tubes, Erlenmeyer flasks, or other containers for the purpose of controlling their temperature.

We do not know which substances are contained within these vessels.

Many substances are:

- **flammable, easily ignited or explosive**
- **hazardous to health**
- **environmentally unsafe**

i.e.: **dangerous**

The user alone is responsible for the handling of these substances!

- Always properly seal all sample containers.



Notice:

There is a danger of electrochemical oxidation or corrosion when using test-tube racks or samples made of non-ferrous metal.

- Avoid using these types of racks or samples.
- Use only original JULABO test-tube racks.

6.2. Bath fluids

Recommended bath fluids: soft/decalcified water.



Caution:

Poor water quality may result in corrosion in the bath.

The quality of water (tap water) depends on local conditions.

- Due to the high concentration of lime, hard water is not suitable for temperature control because it leads to calcification in the bath.
- Ferrous water can cause corrosion - even on stainless steel.
- Chloric water can cause pitting corrosion.
- Distilled and deionized water is unsuitable. Their special properties cause corrosion in the bath, even in stainless steel.

JULABO takes no responsibility for damages caused by the selection of an unsuitable bath fluid.

Please contact JULABO before using other than recommended bath fluids.

Do not use flammable bath fluids!

6.3. Filling / Draining

Filling

- Recommended maximum filling level is 25 mm below the tank rim
- Minimum filling level: approx. 1 cm above the perforated stainless steel base plate.




Note:

The working filling level depends on size and number of the items (fixtures) to be placed inside.

The recommended procedure is to fill the water bath only partially, place the items (fixtures) inside and then correct the filling level (adding or removing liquid) as required.

Draining

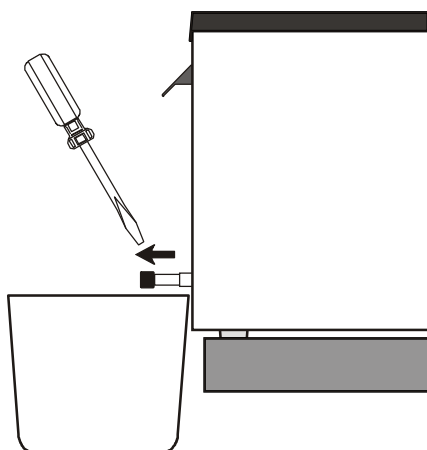

- Press the mains switch to turn the water bath off
- TW2
Take the water bath in both hands and pour out the bath liquid into a suitable vessel.
- TW8, TW12, TW20
Place a suitable collecting bucket or tub underneath the unit for draining the used bath liquid.
- To drain the liquid open the drainage screw (11) on the side of the water bath.
- After the liquid has been fully drained, securely tighten the drainage screw (11) again.



Warning:

There are thermal dangers when opening the bath cover:
Burn, scald, hot steam, hot parts and surfaces that can be touched.

- Do not drain the bath fluid while it is hot!
- Check the temperature of the bath fluid prior to draining (by switching the unit on for a short moment, for example).

Recommendation:

Use the water bath cover to keep temperature losses to a minimum. This is especially important for working temperatures above 60 °C.

6.4. Maintaining a constant water level / Counter cooling

For cooling tasks near the ambient air temperature the liquid level/cooling set can be used for counter cooling.

By special pipe routing, cool faucet water is continuously supplied to the water bath, while at the same time, the heated water is drained via the overflow connection of the Level/Cooling set.

ⓘ A specific water flow rate of 100 ml/minute is sufficient to compensate for the characteristic temperature.



Caution:


Securely attach all tubing to prevent slipping. Observe the laws and regulations of the water distribution company valid in the location where the unit is operated.

Use of the liquid level/cooling set for a continuous supply of faucet water:

1. to keep the water level constant, especially for applications up to the boiling point (supply of faucet water only in the amount of evaporation losses).
2. for counter cooling of cooling tasks near the ambient surrounding temperature (cool faucet water is continuously supplied to the water bath, while at the same time, the heated water is drained via the overflow connection of the liquid level/cooling set).

Liquid level/cooling set


Order number: 8 970 415

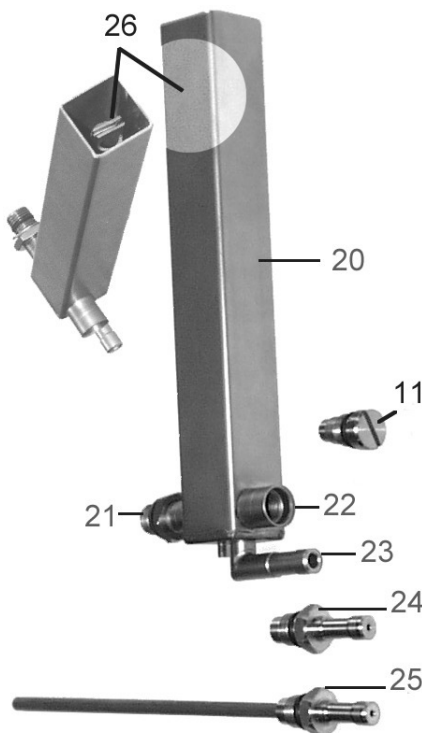

- 11 drainage screw on water bath
- 20 compensation reservoir
- 21 connecting sleeve
- 22 supply/drainage sleeve
- 23 overflow sleeve
- 24 adaptor screw for constant liquid level function
- 25 adaptor screw assy. for counter cooling function and simultaneous constant liquid level control
- 26 adjuster screw for filling level adjustment



6.5. Accessoires

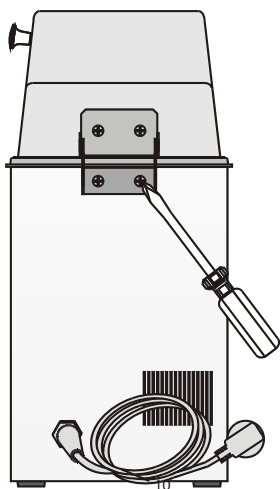

Lift-up Makrolon[®] covers

Order numbers:

TW2	TW8	TW12	TW20
8 970 289	8 970 286	8 970 287	8 970 288

The lift-up Makrolon[®] covers are supplied with pre-assembled hinges. Screw the hinges with the included screws to the rear side of the bath by means of a screwdriver.

TW2

Insert capacity for test tube racks to 100 °C, stainless steel

TW2 1

Order No.	Test tube rack	
8 970 330	for 24 test tubes	16/17 mm dia.

TW8

TW12

TW20

Insert capacity for test tube racks to 100 °C, stainless steel

TW8	TW12	TW20
2	3	4

Order No.	Test tube racks	
8 970 344	for 50 test tubes	16/17 mm dia.
8 970 345	for 90 test tubes	12/13 mm dia.
8 970 346	for 90 microliter tubes	16/17 mm dia.
8 970 347	for 21 test tubes	30 mm dia.

TW8

TW12

TW20

Insert capacity for test tube racks to 80 °C, Polypropylene

TW8	TW12	TW20
2	3	4

Order No.	Test tube racks	
8 970 380	for 60 test tubes	16/17 mm dia.
8 970 381	for 90 test tubes	12/13 mm dia.
8 970 382	for 90 microliter tubes	16/17 mm dia.
8 970 383	for 21 test tubes	30 mm dia.

Cooling installation / continuo water supply

Recommendation: for model TW8, TW12, TW20
For continuous water supply and counter-cooling

Order No.	Description
8 970 415	Liquid level/cooling set
8 970 416	Cooling coil

Water bath protective media

We recommend the use of the "Aqua-Stabil" protective media to eliminate the formation of algae, bacteria, and other micro-organisms.

Order No.	Description
8 940 006	6 bottles 100 ml each
8 940 012	12 bottles 100 ml each

7. Operating procedures

7.1. Power connection



Caution:

- ConnOnly connect the unit to a power socket with earthing contact (PE – protective earth)!
- The power supply plug serves as safe disconnecting device from the line and must be always easily accessible.
- Never operate equipment with damaged mains power cables.
- Regularly check the mains power cable for material defects.
- We disclaim all liability for damage caused by incorrect line voltages!

Check to make sure that the line voltage matches the supply voltage specified on the identification plate.
Deviations of $\pm 10\%$ are permissible.

7.2. Switching on / Start - Stop




- **Switching on:**

Turn on the mains power switch (1).

- ① The unit performs a self-test. All segments of the 4-digit LED temperature DISPLAY and all indicator lights will illuminate. Then the software version (example: n 11.0) appears.

Together with the display of the water bath temperature the operating state is also displayed. (Example: 18.5 °C)






The heat-up phase is indicated by the yellow heating control light  which will blink in regular intervals when the selected temperature has been attained.

- **Switching off:**

Turn the unit off with the mains power switch.

7.3. Setting the temperatures

- ① Setting can be carried out in the start/stop condition.

1. Press one of the keys   for a short moment. The setpoint value instead of the actual value is indicated on the display for about 8 seconds. The value can now be changed.
2. Change value:
Press  to set a higher value.
Press  to set a lower value.
Keep the keys depressed for the value to change fast.
3. Press enter  to store the value.




Notice:

When the working temperature is higher than 50 °C, it might happen that due to strong production of steam there is considerable dripping on the inside of the lift-up Makrolon[®] cover. Some drops may fall directly into the material to be tempered.

- Always properly seal all sample containers.

8. Troubleshooting guide / Error messages

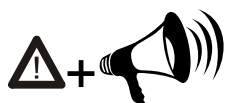


Whenever the microprocessor electronics registers a failure, a complete shutdown of the heater and circulating pump is performed. The alarm light "⚠" illuminates and a continuous signal tone sounds. The LED temperature display indicates the cause for the alarm in form of a code.



- Safety sensor or working temperature sensor.
- The water bath is operated without bath fluid, or the liquid level is insufficient.
Replenish the bath tank with the bath fluid.



- Cable of the working temperature sensor interrupted or short-circuited.

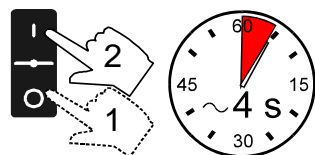

After eliminating the malfunction, press the mains power switch off and on again to cancel the alarm state. If the unit cannot be returned to operation, contact an authorized JULABO service station.

9. Cleaning / repairing the unit



Caution:

Improper maintenance or repair can result in electric shock or damage to the unit.

- Repairs and any other work are to be carried out only by qualified service personnel authorized by JULABO LABORTECHNIK GmbH.
- Always turn off the unit and disconnect the mains cable from the power source before performing any service or maintenance procedures, or before moving the unit.
- Prevent humidity from entering into the water bath.
- Do not use alcohol-based or solvent-based cleaning agents. These cleaning agents will result in damage and cracks in the Makrolon® cover.

- Cleaning:** For cleaning the bath tank and the immersed parts of the water bath, use low surface tension water (e.g., soap suds). Clean the outside of the unit using a wet cloth and low surface tension water.
- The JULABO Water Baths are designed for continuous operation under normal conditions. Periodic maintenance is not required.
- The tank should be filled only with a bath fluid recommended by JULABO. To avoid contamination, it is essential to change the bath fluid from time to time.
- Repairs:** Before asking for a service technician or returning a JULABO instrument for repair, please contact an authorized JULABO service station.
- Returning a unit:** When returning the unit:
- Clean the unit and, if necessary, decontaminate the unit in order to avoid endangering service personnel.
 - Attach a short fault description.
If you intend to return your JULABO unit to us, you will find a Service Return Form on our website www.julabo.de. Please use this as a delivery note and include it to the unit or send it in advance either by Fax or E-Mail.
 - During transport the unit has to stand upright. Mark the packing correspondingly.
 - When returning a unit, take care of careful and adequate packing.
 - JULABO is not responsible for damages that might occur from insufficient packing.



JULABO reserves the right to carry out technical modifications with repairs for providing improved performance of a unit.